

# PRODUCT CATALOGUE

## PRIORITY PROTECT MASK



A black and white photograph of a man with curly hair wearing a dark, textured, reusable protective mask. The mask has a small rectangular detail on the side. The background is a plain, light color.

## **PRODUCT BRIEF**

There are different masks available in the market now. However, it is important to be able to differentiate between a protective mask for you and a surgical mask meant for the doctors.

What we need is a reusable protective mask with an anti-virus, antibacterial, anti-dust, anti-pollution mask with 6 layer particle filtration system. A Priority Protect mask is engineered exactly with these features, also keeping in mind comfort, breath-ability and safety. Priority's Protect Mask is engineered with a 6 Layer Filtration System. First filtration Layer helps to filter coarser unwanted particles, middle filtration is of spun and melt blown which works as bacteria & small particles filter & prevents any bacteria sized 0.3 micron, last filtration layer is a cotton fabric for moisture control and has an anti-bacterial feature.

# PRODUCT FEATURES

Minimal  
eyewear  
fogging


Adjustable  
earloops

Contour  
fitting design

Flawless  
construction

Wide face  
coverage  
for maximum  
protection

Engineered  
6 layered  
particle  
filtration

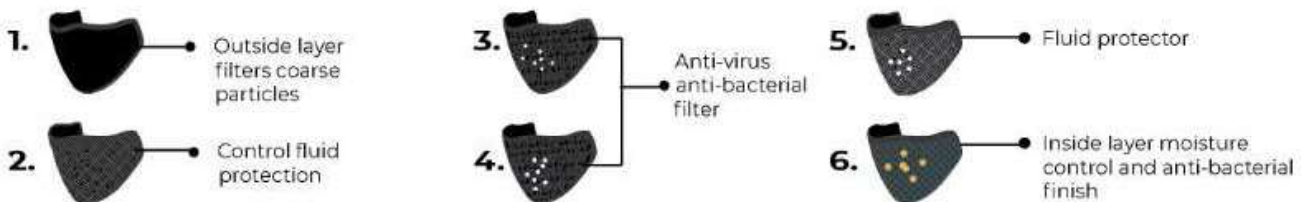

## Why Priority Protect Mask

- Easy to store when not being used
- Fits every face shape and size with its 2 panel design
- Crafted with a six layer particle filtration system
- Equipped with soft elasticated ear loops
- Discards exhaled air and heat build up
- Allows easy breathing and facial moments

# **PRODUCT TESTING AND CERTIFICATION**

Priority Protect Mask has been tested & certified as per standards by ISO 9001:2015, ISO 16603:2004, ISO 13485:2016, CE COMPLIANCE, GMP COMPLIANCE. Our protect mask delivers exceptionally high level of Protection. The test details are below:


1. 95%+ bacterial filtration efficiency (BFE) for over 3um (micron) sized bacteria.
2. 95%+ Particle filtration efficiency (PFE) for over 3um (micron) sized bacteria.
3. Splash Resistance (Hydrostatic Test at 160 mmHg) which protects against accidental fluid/sneeze droplets from reaching user's nose or mouth.
4. Easy to breathe & wear (Differential Pressure test {Delta P} with <35 pa/cm square)



# PRODUCT USAGE

Designed with a six layer filtration system, Priority Protect mask is conducive for respiratory protective the city air infectious and harmful particles like dust, sand, pollen, smog and other infections.

This reusable and washable mask that is cost-effective, is a definite way to protect your lungs.

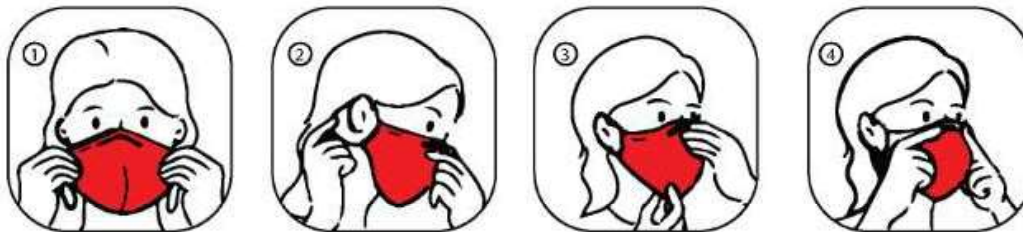

1. Hold the ear strings of the mask and stretch it out, making sure the concave side is facing your mouth and that the top line comes above your nostrils.
2. Wear the mask, ensuring that it covers your nose, mouth and chin, and secure the ear strings around your ears.
3. Adjust the earloops strap for the right and comfortable fit and ensure your safety.

## **Note:**

1. In case of difficulty in breathing, please change the mask.
2. Store the Priority Protect mask away from the sunlight, at a dry place.
3. Wear the mask properly as per the instructions or it may lead to poor performance of the mask and your protection may be hampered.